


2-Oxiranecarboxylicacid, 3-propyl-, methyl ester (+/-) Chiral Separation on MEGA-DEX Columns

Column: **MEGA-DEX DMT Beta - 0.25mm, 0.25 μ m, 25m**

Catalog Code: DEX-DMT-B-025-025-25

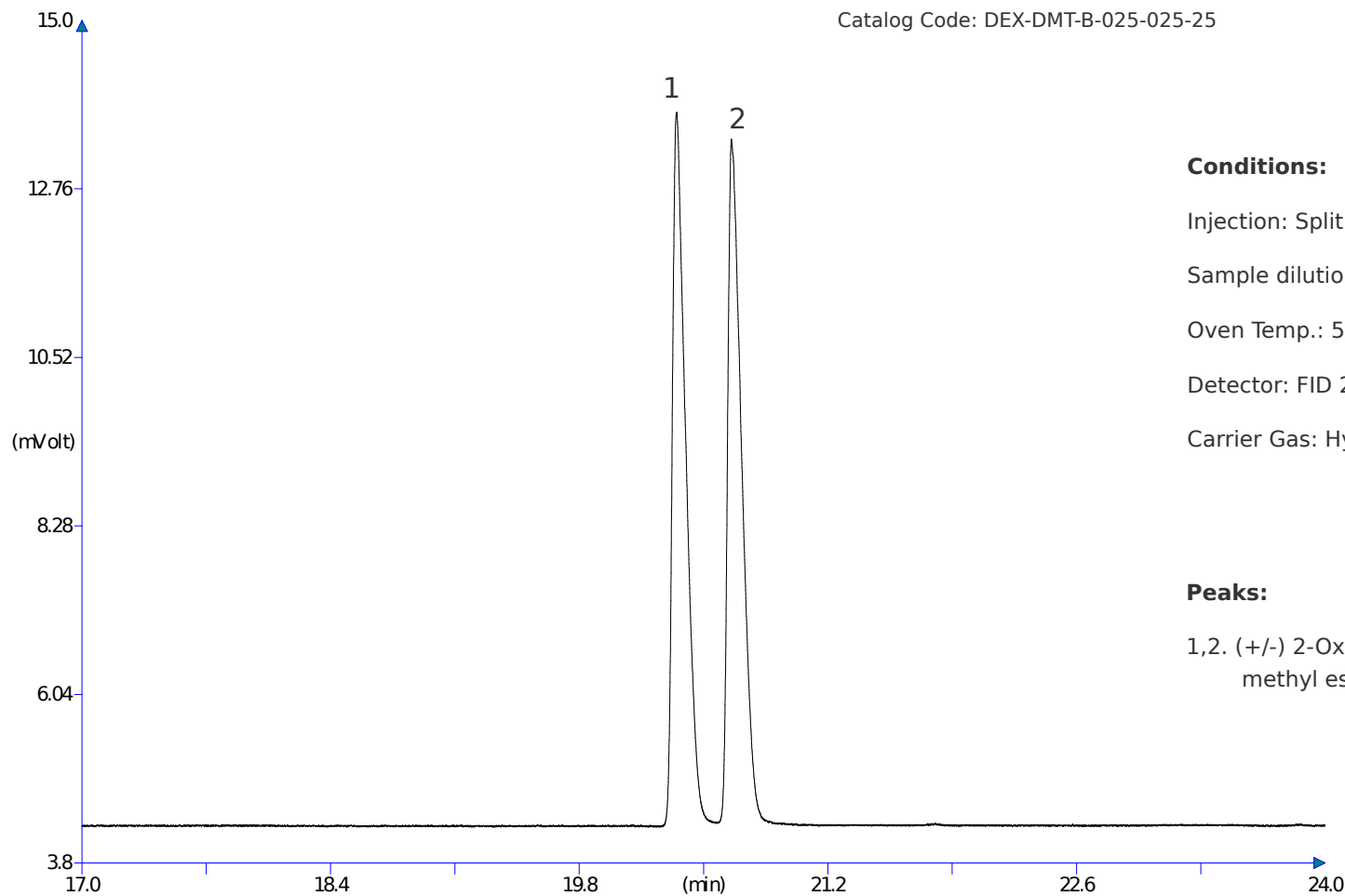

Conditions:

Injection: Split 250°C, 0.5 μ L, 90mL/min Split Flow.

Sample dilution: 1mg/1mL in Methylene Chloride.

Oven Temp.: 50°C, 2°C/min, 200°C.

Detector: FID 250°C.

Carrier Gas: Hydrogen, 70 kPa.

Peaks:

1,2. (+/-) 2-Oxiranecarboxylicacid, 3-propyl-,
methyl ester

