


Acetylated *trans*-Cyclohexanedimethanol (+/-) Chiral Separation on MEGA-DEX Columns

Column: **MEGA-DEX DMT Beta - 0.25mm, 0.25 μ m, 25m**

Catalog Code: DEX-DMT-B-025-025-25

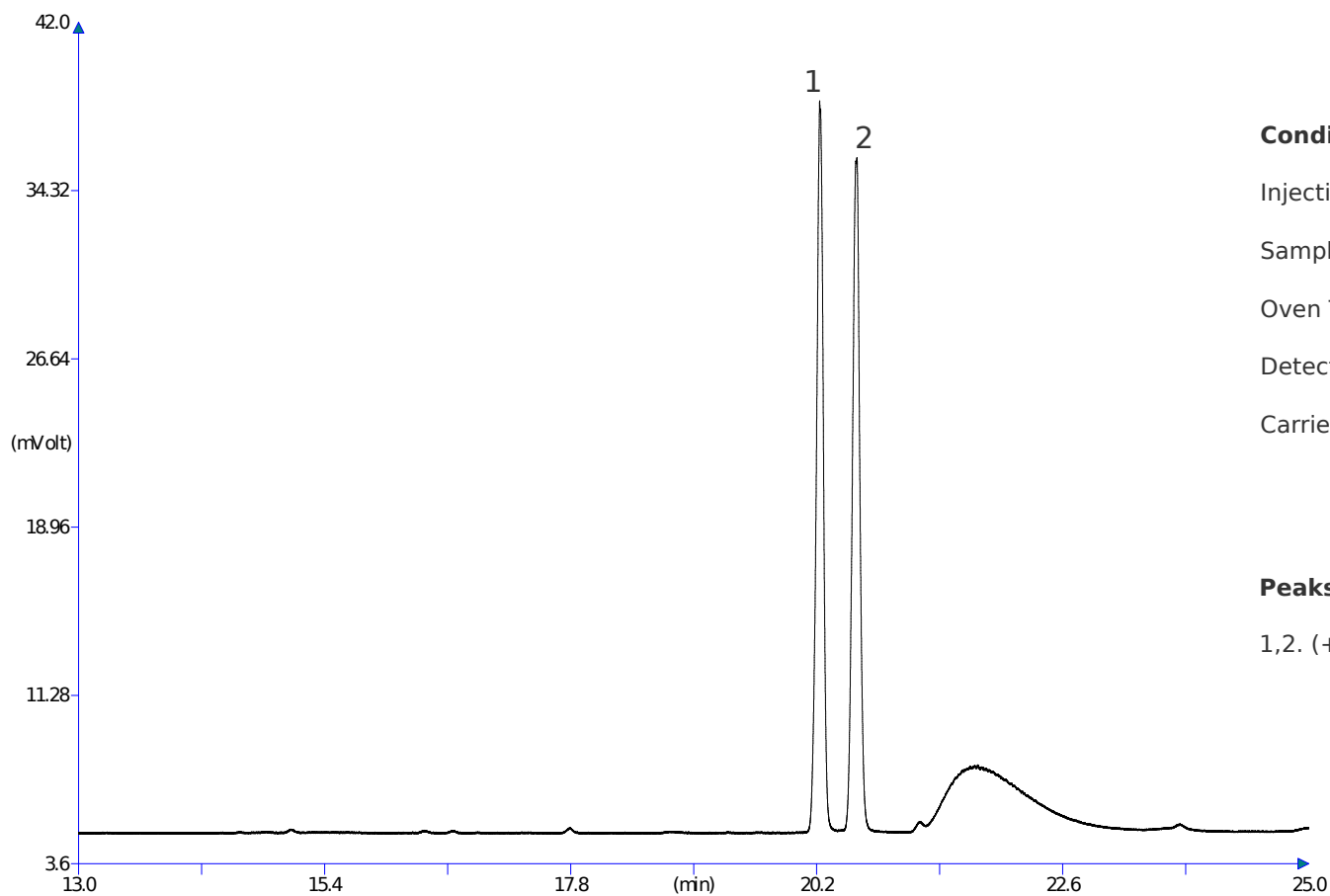

Conditions:

Injection: Split 250°C, 0.5 μ L, 100mL/min Split Flow.

Sample dilution: 1mg/1mL.

Oven Temp.: 120°C, 2.5°C/min, 200°C.

Detector: FID 250°C.

Carrier Gas: Helium, 110 kPa.

Peaks:

1,2. (+/-) *trans* Cyclohexanedimethanol

