


Alcohol Ethoxylates

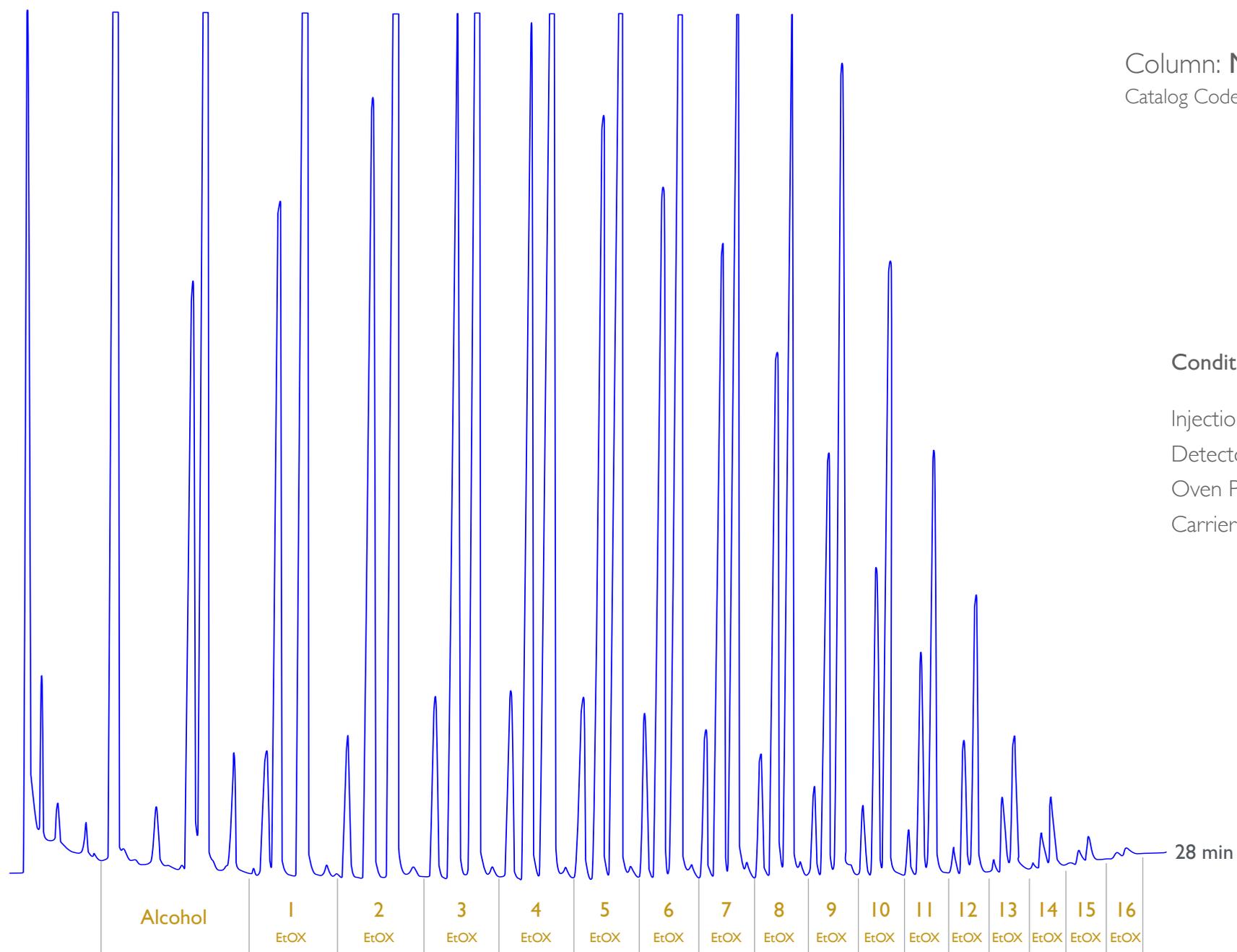

Column: **MEGA-SE52 HT** - 0.32mm, 0.10 μ m, 10m
Catalog Code: HT-SE52-032-010-10

Conditions:

Injection: cold On-Column,

Detector: FID 390°C.

Oven Program: 65°C (1min), 10°C/min up to 380°C.

Carrier Gas: Helium, 2mL/min.



since
1980



MEGA[®]
improve your GC analysis

contact us: info@mega.mi.it



follow us [@MEGAColumns](https://twitter.com/MEGAColumns)

www.mega.mi.it